

Respimat[®] Soft Mist Inhaler


Use it correctly for better medicine delivery and better breathing!



LOADING THE MEDICATION

1 Remove cap

Press the safety catch and pull off the clear base.

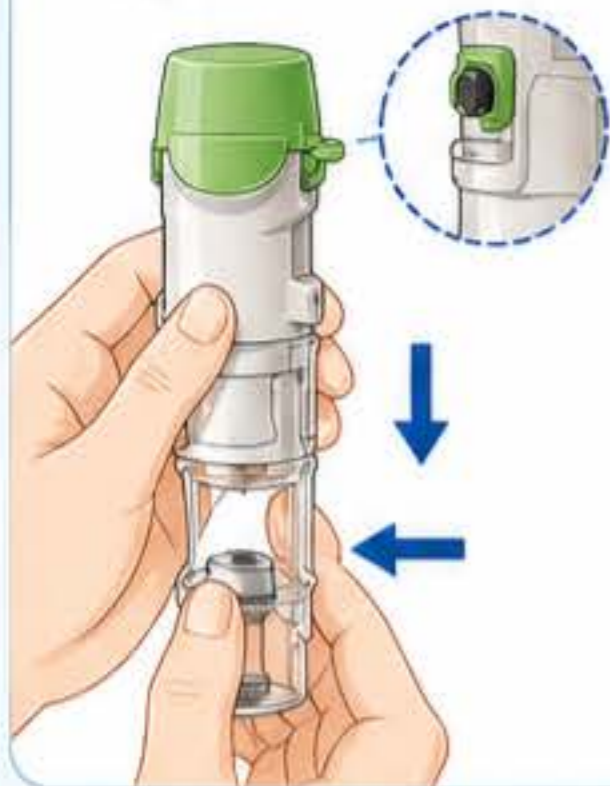

2 Insert cartridge

Place the narrow end of the cartridge into the inhaler.



3 Push and click

Push the cartridge firmly into the inhaler until it clicks. Once inserted, **do not remove it.**



4 Prime

Reattach the clear base. After first loading, press the dose-release button **3 times** to prime until you see a mist before first use.



USING YOUR INHALER

1 Rotate

Hold the inhaler upright with the green cap closed. Turn the clear base in the direction of the red arrows on the label until it clicks.



2 Open

Open the green cap until it snaps fully open.



3 Breathe out

Breathe out slowly and fully, away from the inhaler.



4 Inhale & press

Place the mouthpiece between your lips, without blocking the air vents. Point the inhaler toward the back of your throat. Inhale slowly and deeply through your mouth while pressing the dose-release button. Continue to inhale slowly and deeply. **Hold your breath for 10 seconds.**





5 Close


Close the green cap until it snaps back into place.



TIPS FOR BEST RESULTS

 Close your lips around the mouthpiece but do not block the air vents.

 Inhale slowly and deeply to allow the medicine to reach deep into your lungs.


 Use your inhaler at the same time every day. Close the green cap and keep it clean.


IMPORTANT NOTES

- Once the cartridge is inserted into the inhaler, do not remove it.
- If not used for 3 days or more, prime once (spray into the air **1 time**).
- Clean the mouthpiece regularly: Wipe the mouthpiece with a damp cloth or tissue (do not remove). Let it air dry completely to avoid blocking the mist.



STORAGE

 Store at room temperature. Keep away from direct sunlight and high heat.

 Keep out of reach of children.